

CANTY

PROCESS TECHNOLOGY

MACROFLOW™ PARTICLE SIZING SYSTEM



THE CANTY ADVANTAGE

The Canty Macroflow™ is a unique vision system designed to analyze particles in suspension on a macro scale using CantyVision™ software. Particles are analyzed in real time for two dimensional size, shape, and color. A sample or continuous flow is analyzed under actual operating conditions of temperature and pressure. The image sensor is a high resolution, high speed CCD camera coupled to an adjustable macro-zoom lens system to cover a wide range of particle size. Ethernet controlled zoom and focus capabilities with preset calibration options for multiple lens settings are offered. Canty Fuseview technology is incorporated into the structure. Fused glass far exceeds all other windows in safety and performance. Canty fiberoptic cold lighting is integrated into the Macroflow™ design, so no additional heat is introduced into the process. The rugged construction allows for use in the field as well as in the lab. The flush design has no pockets and provides a seamless flow pattern where particles can be analyzed without concerns of attrition or damaging effects.

APPLICATIONS

- Measure dissolving rates for soluble material
- Polymer Bead Growth and Shape Analysis
- Mining and Aggregate Particle Size
- Color Analysis for Lab or In-line
- In-Line Particle Size for Process Control
- Gel, Emulsion and Bubble Classifier
- Micro-Reactor Analysis - Dual Inputs and Flush Loop with Sanitary Connections
- Percent Solids / Turbidity
- End Product Quality Control and Lot Certification
- Replace Laser and Sieve with Two Dimensional Distribution by Shape

FEATURES

- Particle Size, Shape and Distribution with Real Time Analysis
- Fused Glass Pressure Barrier Standard
- Rated 150 PSI @ 500 Deg. F., Options Through 6,000 PSI.
- On-line or Lab Plus (remote control)
- Seamless Design. Unit Can Be Orientated Horizontal or Vertical
- In-Line Dilution Port
- Instantaneous Mixing Port
- Regulated Temperature Option
- Bypass Loop Option Available for Flow Control
- 3/4 Inch Sanitary Connections Standard

CANTY JM Canty Inc | Buffalo, NY USA
JM Canty Intl Ltd | Dublin, Ireland

Ph: (716) 625 4227
Ph: + 353 (01) 882 9621

Fax: (716) 625 4228
Fax: +353 (01) 882 9622

www.jmcanty.com

P/N TA10491-100 Rev. 1

Ordering Information

HOW TO ORDER: Select the appropriate symbols and build a model number as shown:

EXAMPLE:

MAC - L A 1 K 1 1 B 2 1 A Z

TYPE _____

- L - Laminar Flow
- H- High Pressure

SYSTEM VIDEO FORMAT _____

- A - EIA (US B&W) 640 x 480
- B - NTSC (US Color) 640 x 480
- C - CCIR (European B&W) 768 x 576
- D - PAL (European Color) 768 x 576
- E - Ethernet (Color) 1600 x 1200, 800 x 600
- F - Ethernet (B&W) 1600 x 1200, 800 x 600
- G - Gigabit Ethernet (B&W) 1600 x 1200
- H - Gigabit Ethernet (Color) 1600 x 1200

ENVIRONMENTAL RATING _____

- 1 - Lab unit (General purpose enclosure)
- 2 - NEMA 4 weatherproof, IP 66
- 3 - Lab Plus (Ethernet remote control)
- 4 - Explosion proof (US-120 V)
or Flameproof (Europe-230 V)

SYSTEM ILLUMINATION _____

- K - Backlight
- L - Darkfield
- M - Cross Polarized
- N - Amplified / Focused
- P - Phase Contrast

PARTICLE SIZE RANGE _____


P.S. View	P.S. Measure	
1 - (4.8µm - 6,000µm)	(10µm - 1,200µm) w/ 1/4" Deep flowpath @ 1600 x 1200	
2 - (4.8µm - 12,000µm)	(10µm - 2,400µm) w/ 1/2" Deep flowpath @ 1600 x 1200	
3 - (10µm - 6,000µm)	(20µm - 1,200µm) w/ 1/4" Deep flowpath @ 800 x 600	
4 - (10µm - 12,000µm)	(20µm - 2,400µm) w/ 1/2" Deep flowpath @ 800 x 600	
5 - (13µm - 6,000µm)	(25µm - 1,200µm) w/ 1/4" Deep flowpath @ 640 x 480, 768 x 576	
6 - (13µm - 12,000µm)	(25µm - 2,400µm) w/ 1/2" Deep flowpath @ 640 x 480, 768 x 576	

For smaller sizes refer to Microflow sheet TA9309-1.

All registered trademarks are the property of their respective owners.



Flow Diagram

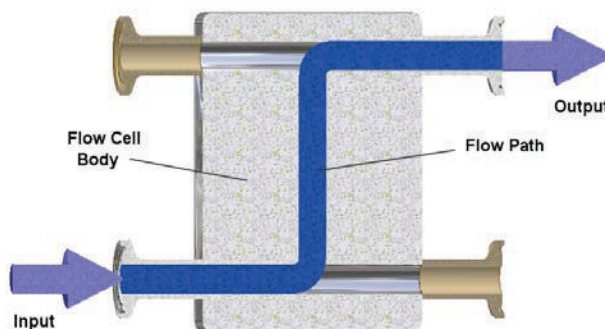

SPECIFICATIONS

- Power: 120V AC / 60 Hz
(230V AC / 50Hz, models available)
- 3/4 inch Sanitary flow connections (two in and two out, with two plugs)
- Weatherproof Seals on Camera and Light Source
- Optional Laboratory Flow Loop Sample Feed
- Gasket Materials: Teflon®

<p>CANTY JM Canty Inc JM Canty Intl Ltd</p>	<p>Buffalo, NY USA Dublin, Ireland</p>	<p>Ph: (716) 625 4227 Ph: + 353 (01) 882 9621</p>	<p>Fax: (716) 625 4228 Fax: +353 (01) 882 9622</p>
--	--	---	--